



Anti 4R-tau

BACKGROUND

Tau is a microtubule-associated protein. In adult human brain, 6 tau isoforms are expressed, which are divided into two groups, 3-repeat (3R) and 4-repeat (4R) tau isoforms. Tau assembles into insoluble, filamentous, and hyperphosphorylated inclusions in many neurodegenerative diseases, and different tau isoforms are accumulated in the diseased brains, namely, 6 tau isoforms in AD, 3R tau isoforms in Pick's disease, and 4R tau isoforms in PSP and CBD. Therefore, isoform-specific tau antibodies are useful tools for immunohistochemical and biochemical studies of tau species in diseased brains. Tau deposited in AD is extensively deamidated at N279. This anti-4R specifically recognizes 4R tau isoforms regardless of deamidation and strongly stained tau in AD brain.

Product type	Primary antibody
Immunogen	synthetic N279D 4R-tau peptide
Raised in	(VQIIDKKLDLSNVQSKC) Rabbit
Myeloma Clone number	-
Isotype	-
Source	-
Purification	Anti-serum
Form	-
Concentration	Liquid. Anti-serum with 0.1% NaN3 as a preservative
Volume	-
Label	50 uL
Specificity	Unlabeled
Cross reactivity	4R-tau isoforms
Storage	Human, Mouse, Rat Store below -20°C. (below -70°C for prolonged storage). Aliquot to avoid cycles of freeze/thaw.


Application notes ▪ Western blotting: 1/2000 - 1/5000 (Fig. 1)
Recommended dilutions ▪ Immunohistochemistry: 1:3000 - 1/20000 (Fig. 2)
 a higher dilution is recommended for immunohistochemistry.
 Other applications have not been tested.
 Optimal dilutions/concentrations should be determined by the end user.

References

- 1) Hasegawa M, Watanabe S, Kondo H, et al, 3R and 4R tau isoforms in paired helical filaments in Alzheimer's disease. Acta Neuropathol. 2013 Nov 9
- 2) Dan A, Takahashi M, Masuda-Suzukake M, et al. Extensive deamidation at asparagine residue 279 accounts for weak immunoreactivity of tau with RD4 antibody in Alzheimer's disease brain. Acta Neuropathol



ANTIBODY CHARACTERIZATION

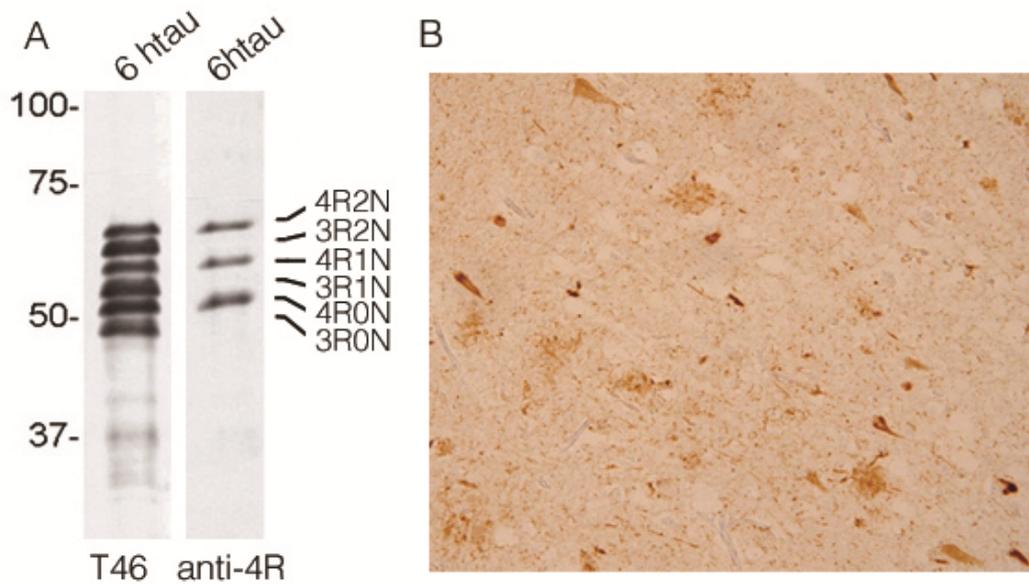

(Fig1)

Fig1: Immuno blot analysis for 6 kinds of tau isoforms (3-repeat and 4-repeat) with anti total tau (T46) and Anti 4R (TIP-4RT-P01).

Fig2: Immuno staining for paraffin embedded section of cerebral cortex prepared from AD brain (retrieved with proteinase K and formic acid)

(Fig2)



COSMO BIO Co., LTD.

【JAPAN】
TOYO EKIMAE BLDG. 2-20, TOYO 2-CHOME,
KOTO-KU. TOKYO 135-0016, JAPAN
Phone: +81-3-5632-9610
FAX: +81-3-5632-9619
URL: <https://www.cosmobio.co.jp/>



COSMO BIO USA

【Outside Japan】
2792 Loker Ave West, Suite 101
Carlsbad, CA 92010, USA
email: support@cosmobiousa.com
URL: www.cosmobiousa.com
Phone/FAX: (+1) 760-431-4600