



MONOCLONAL ANTIBODY

For research use only. Not for clinical diagnosis.

Catalog No. NM-MA-004

Anti- 5-Fluorouracil

purified with protein G affinity column


BACKGROUND

5-Fluorouracil (5-FU) is a pyrimidine analogue and inhibits an enzyme called thymidylate synthetase, which results in inhibition of DNA replication. Thus, 5-FU is used as a drug in the treatment of cancers including colorectal cancer, pancreatic cancer and skin cancer.

Product type	Primary antibody
Immunogen	5-Fluorouridine—BSA (5-FU-BSA) [5-FU]/[BSA]=10.2
Myeloma	P3-X63-Ag8.653
Clone number	H3-17
Isotype	IgG1, λ
Host	Mouse
Purification	Protein G affinity column
Form	Liquid
Storage buffer	PBS, No preservative is contained.
Concentration	100 μ g/mL
Volume	50 μ L
Label	Unlabeled
Specificity	5-FU (both free 5-FU and protein-bound 5-FU)
Cross reactivity	The antibody does not have cross-reactivity with Uracil, Cytosine, Thymine, 5,6-Dihydro-5-methyluracil and 5,6-Dihydrouracil.
Storage	Lyophilized form: store at -20 to -80 °C. Reconstituted form: store at -20 °C. After reconstitution, it is stable for at least 1 year when stored at -20 °C. It should be divided into small quantity to avoid many freezing and thawing.

Other

< 5-Fluorouracil >

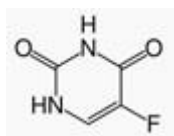

CAS.No.	51-21-8
Molecular Formula	C ₄ H ₃ FN ₂ O ₂
Molecular Weight	130.08

Application notes	ELISA
Recommended dilutions	• ELISA : 1:100-1000

Other applications have not been tested.

Optimal dilutions/concentrations should be determined by the end user.

ANTIBODY CHARACTERIZATION

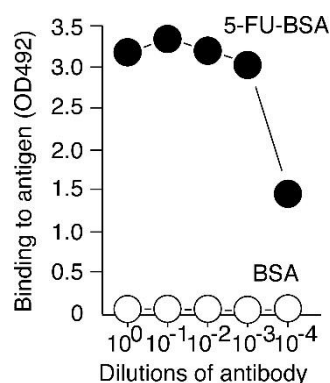


Fig.1 Monoclonal antibody (H3-17) shows high binding to 5-FU-BSA but undetectable binding to BSA. Different dilutions of antibody were tested for binding to immobilized antigens (100 ng / well) in a direct ELISA.

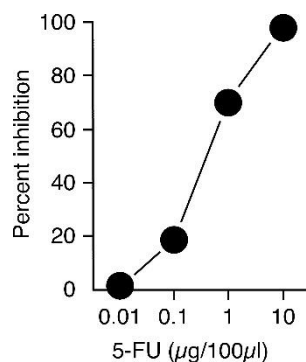


Fig.2 Monoclonal antibody (H3-17) is capable of binding to free 5-FU.. Free 5-FU efficiently inhibits the antibody binding to immobilized 5-FU-BSA (5 ng / well), which was detected by a competitive ELISA.

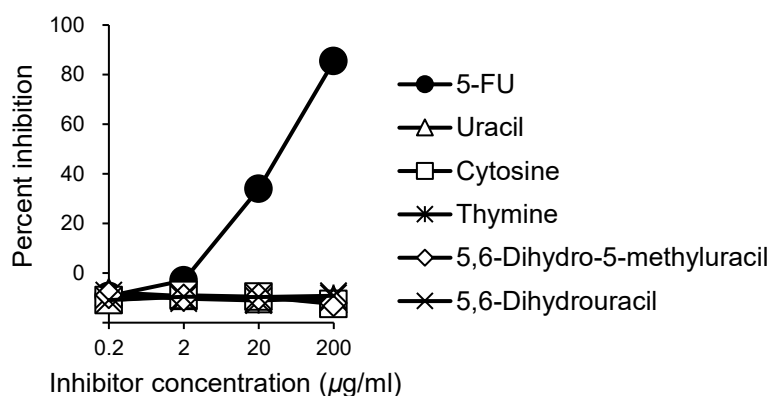


Fig.3 Binding specificity of Anti-5-FU antibody. Anti-5-FU monoclonal antibody (H3-17) binds to 5-FU, but not to Uracil, Cytosine, Thymine, 5,6-Dihydro-5-methyluracil and 5,6-Dihydrouracil.

RELATED PRODUCTS

Product Name	Maker	Cat#
5-FU-BSA (5-Fluorouracil Bovine Serum Albumin conjugate)	CSR	NM-MA-R001

ELISA Protocols

A. Direct ELISA (Fig. 1)

- 1) Prepare 5-FU-BSA solutions in PBS at the concentration of 2 µg / mL.
- 2) Distribute 50 µL / well of the 5-FU-BSA solution to 96 well microtiter plates.
- 3) Seal the plates with plate seals, and leave overnight at 4 °C.
- 4) Wash the 5-FU-BSA-coated plates 5 times with 150 µL / well of PBS-T (0.05% Tween-20 in PBS).
- 5) Distribute 150 µL / well of 2% FBS in PBS to each well to prevent non-specific antibody binding.
- 6) Incubate 30 minutes at 37 °C.
- 7) Wash the plates 5 times with 150 µL / well of PBS-T.
Discard PBS-T, and tap the plate on the paper towel to remove the solution completely.
- 8) Prepare serial dilutions of H3-17 antibody solutions in PBS.
- 9) Distribute 100 µL / well of H3-17 antibodies and incubate 30 minutes at 37 °C.
- 10) Wash the plates 5 times with 150 µL / well of PBS-T.
Discard PBS-T, and tap the plate on the paper towel to remove the solution completely.
- 11) Distribute 100 µL / well of 1:2500 HRP-goat anti-mouse IgG (H+L) (invitrogen, Cat. No. 62-6520) diluted with PBS to each well and incubate 30 minutes at 37 °C.
- 12) Wash the plates 5 times with 150 µL / well of PBS-T.
Discard PBS-T, and tap the plate on the paper towel to remove the solution completely.
- 13) Distribute 100 µL / well of the substrate solution [o-Phenylene diamine 8 mg, H₂O₂ (30%) 4 µL, Citrate-phosphate buffer (pH5.0) 20 mL] to each well and incubate 30 minutes at 37 °C.
- 14) Add 50 µL / well of 2M H₂SO₄ to each well and stop enzyme reaction.
- 15) After gentle mixing, determine the absorbance at 492 nm of each well by a spectrophotometer.

B. Competitive ELISA (Fig. 2)

- 1) Prepare 5-FU-BSA solutions in PBS at the concentration of 0.1 µg / mL (5 ng / well).
- 2) Distribute 50 µL / well of the 5-FU-BSA solution to 96 well microtiter plates.
- 3) Seal the plates with plate seals, and leave overnight at 4 °C.
- 4) Wash the 5-FU-BSA-coated plates 5 times with 150 µL / well of PBS-T (0.05% Tween-20 in PBS).
Discard PBS-T, and tap the plate on the paper towel to remove the solution completely.
- 5) Distribute 150 µL / well of 2% FBS in PBS to each well to prevent non-specific antibody binding.
- 6) Incubate 30 minutes at 37 °C.
- 7) Wash the plates 5 times with 150 µL / well of PBS-T.
Discard PBS-T, and tap the plate on the paper towel to remove the solution completely.
- 8) Prepare 5-FU (competitor, 50 µL) solutions in tubes which concentrations are 0, 0.01, 0.1, 1, 10 ug/ 50 µL PBS.
Add 50 µL of 1:500 H3-17 antibody solution to each tube, which gives 50% of the maximum binding to the solid-phase antigen. And mix gently.
- 9) Distribute 100 µL /well of mixtures to each well and incubate 30 minutes at 37 °C.
- 10) Wash the plates 5 times with 150 µL / well of PBS-T.
Discard PBS-T, and tap the plate on the paper towel to remove the solution completely.
- 11) Distribute 100 µL / well of 1:2500 HRP-goat anti-mouse IgG (H+L) (invitrogen, Cat. No. 62-6520) diluted with PBS to each well and incubate 30 minutes at 37 °C.
- 12) Wash the plates 5 times with 150 µL / well of PBS-T.
Discard PBS-T, and tap the plate on the paper towel to remove the solution completely.
- 13) Distribute 100 µL / well of the substrate solution [o-Phenylene diamine 8 mg, H₂O₂ (30%) 4 µL, Citrate-phosphate buffer (pH5.0) 20 mL] to each well and incubate 30 minutes at 37 °C.
- 14) Add 50 µL / well of 2M H₂SO₄ to each well and stop enzyme reaction.
- 15) After gentle mixing, determine the absorbance at 492 nm of each well by a spectrophotometer.

For research use only, Not for diagnostic use.



COSMO BIO CO., LTD.

【JAPAN】
TOYO EKIMAE BLDG. 2-20, TOYO 2-CHOME,
KOTO-KU. TOKYO 135-0016, JAPAN
Phone: +81-3-5632-9610
FAX: +81-3-5632-9619
URL: <https://www.cosmobio.co.jp/>



COSMO BIO USA

【Outside Japan】
2792 Loker Ave West, Suite 101
Carlsbad, CA 92010, USA
email: support@cosmobiousa.com
Phone/FAX: (+1) 760-431-4600
URL: www.cosmobiousa.com