CircuLex[™]



User's Manual For Research Use Only, Not for use in diagnostic procedures

ELISA Kit for Measuring Human S100A8/MRP8

CircuLex S100A8/MRP8 ELISA Kit

Cat# CY-8061

Intended Use	. 1
Storage	. 1
Introduction	2
Principle of the Assay	2-3
Materials Provided	3
Materials Required but not Provided	4
Precautions and Recommendations	. 5
Sample Collection and Storage	6
Detailed Protocol	. 7-9
Calculations	9
Measurement Range	9
Troubleshooting	9
Reagent Stability	10
Assay Characteristics	10-12
Example of Test Results	
References	

Intended Use

The MBL Research Product **CircuLex S100A8/MRP8 ELISA Kit** is used for the quantitative measurement of human S100A8/MRP8 in serum, plasma, and other biological samples.

Individual users should determine appropriate conditions when using other types of samples.

This assay kit is for research use only and not for use in diagnostic or therapeutic procedures.

Storage

- Upon receipt store all components at 4°C.
- Don't expose reagents to excessive light.



User's Manual



For Research Use Only, Not for use in diagnostic procedures

Introduction

S100A8/MRP8 and its hetero-dimerization partner S100A9/MRP14 belong to the low molecular weight calcium-binding protein, S100 family (1), they are composed of two distinct helix-loop-helix motifs (EF-hands) flanked by hydrophobic regions at either terminus and separated by a central hinge region. In human, 100A8 is usually co-expressed with S100A9. Both proteins are expressed during myeloid differentiation, are abundant in granulocytes and monocytes, and form heterodimeric complexes (2). Although a number of possible functions for S100A8-A9 heterocomplex, including antimicrobial activity, have been proposed, the exact role of these proteins in cell metabolism is still unclear. In human, they have been associated with several inflammatory diseases (3): elevated serum levels of S100A8 and S100A9 have been found in patients suffering from a number of inflammatory disorders including giant cell arteritis (4), cystic fibrosis, rheumatoid arthritis, dermatoses, chronic inflammatory bowel disease, chronic bronchitis (3), some malignancies and autoimmune diseases (5, 6). Both proteins are localized predominantly in the cytoplasm. It was demonstrated with human monocytes that both proteins are secreted by an energy-consuming pathway, which is dependent on an intact microtubule network and involves protein kinase C (7). It was shown that the expression of S100A8 is up regulated by anti-inflammatory mediators (8) and oxidative stress such as UVA and H_2O_2 (9). In addition, this protein seems to be a potent scavenger of hypochlorite (10). Thus this protein may regulate oxidative events associated with inflammation.

Principle of the Assay

The MBL Research Product **CircuLex S100A8/MRP8 ELISA Kit** employs the quantitative sandwich enzyme immunoassay technique. An antibody specific for human S100A8/MRP8 has been pre-coated onto a microplate. Standards and samples are pipetted into the wells and the immobilized antibody binds any human S100A8/MRP8 present. After washing away any unbound substances, an HRP conjugated antibody specific for human S100A8/MRP8 is added to the wells. Following a wash to remove any unbound antibody HRP conjugate, the remaining conjugate is allowed to react with the substrate H_2O_2 -tetramethylbenzidine. The reaction is stopped by addition of acidic solution and absorbance of the resulting yellow product is measured at 450 nm. The absorbance is proportional to the concentration of human S100A8/MRP8. A standard curve is constructed by plotting absorbance values versus human S100A8/MRP8 concentrations of calibrators, and concentrations of unknown samples are determined using this standard curve.



User's Manual



For Research Use Only, Not for use in diagnostic procedures

Summary of Procedure

Add 100 μ L of \perp d	diluted samples to the wells
•	Incubate for 60 minutes at room temp.
Wash ↓ t	he wells
Add 100 μ L of \downarrow H	HRP conjugated anti-human S100A8/MRP8 antibody Incubate for 60 minutes at room temp.
Wash ↓ t	he wells
Add 100 μ L of \downarrow S	Substrate Reagent
Incubate for 1	10-20 minutes at room temp.
Add 100 μ L of \downarrow S	Stop Solution

Measure absorbance at 450 nm

Materials Provided

All samples and standards should be assayed in duplicate. The following components are supplied and are sufficient for the one 96-well microplate kit.

Microplate: One microplate supplied ready to use, with 96 wells (12 strips of 8-wells) in a foil, zip-lock bag with a desiccant pack. Wells are coated with anti-human S100A8/MRP8 antibody as a capture antibody.

10X Wash Buffer: One bottle containing 100 mL of 10X buffer containing Tween[®]-20

Dilution Buffer: One bottle containing 50 mL of 1X buffer; use for sample dilution. Ready to use.

Human S100A8 Standard: One vial, each containing X* ng of lyophilized recombinant human S100A8/MRP8

*The amount is changed depending on lot. See the real "User's Manual" included in the kit box.

20x HRP conjugated Detection Antibody: One vial containing 600 μ L of HRP (horseradish peroxidase) conjugated anti-human S100A8/MRP8 antibody.

Conjugate Dilution Buffer: One bottle containing 12 mL of Conjugate Dilution Buffer.

Substrate Reagent: One bottle containing 20 mL of the chromogenic substrate, tetra-methylbenzidine (TMB). Ready to use.

Stop Solution: One bottle containing 20 mL of 1 N H₂SO₄. Ready to use.





For Research Use Only, Not for use in diagnostic procedures

Materials Required but not Provided

- Pipettors: 2-20 µL, 20-200 µL and 200-1,000 µL precision pipettors with disposable tips.
- Precision repeating pipettor
- Orbital microplate shaker
- Microcentrifuge and tubes for sample preparation.
- Vortex mixer
- (Optional) Microplate washer: Manual washing is possible but not preferable.
- **Plate reader** capable of measuring absorbance in 96-well plates at dual wavelengths of 450 nm/540 nm. Dual wavelengths of 450/550 or 450/595 nm can also be used. The plate can also be read at a single wavelength of 450 nm, which will give a somewhat higher reading.
- (Optional) Software package facilitating data generation and analysis
- 500 or 1,000 mL graduated cylinder.
- Reagent reservoirs
- Deionized water of the highest quality
- Disposable paper towels

CircuLex[™]



For Research Use Only, Not for use in diagnostic procedures

Precautions and Recommendations

- Although we suggest to conduct experiments as outlined below, the optimal experimental conditions will vary depending on the parameters being investigated, and must be determined by the individual user.
- Allow all the components to come to room temperature before use.
- All microplate strips that are not immediately required should be returned to the zip-lock pouch, which must be carefully resealed to avoid moisture absorption.
- Do not use kit components beyond the indicated kit expiration date.
- Use only the microtiter wells provided with the kit.
- Rinse all detergent residue from glassware.
- Use deionized water of the highest quality.
- Do not mix reagents from different kits.
- The buffers and reagents in this kit may contain preservatives or other chemicals. Care should be taken to avoid direct contact with these reagents.
- Do not mouth pipette or ingest any of the reagents.
- Do not smoke, eat, or drink when performing the assay or in areas where samples or reagents are handled.
- Dispose of tetra-methylbenzidine (TMB) containing solutions in compliance with local regulations.
- Avoid contact with the acidic Stop Solution and Substrate Solution, which contains hydrogen peroxide.
- Wear gloves and eye protection when handling immunodiagnostic materials and samples of rat origin, and these reagents. In case of contact with the Stop Solution and the Substrate Solution, wash skin thoroughly with water and seek medical attention, when necessary.
- Biological samples may be contaminated with infectious agents. Do not ingest, expose to open wounds or breathe aerosols. Wear protective gloves and dispose of biological samples properly.
- CAUTION: Sulfuric Acid is a strong acid. Wear disposable gloves and eye protection when handling Stop Solution.

CircuLex[™]



For Research Use Only, Not for use in diagnostic procedures

Sample Collection and Storage

Serum: Allow blood samples to clot for 60 ± 30 minutes. Centrifuge the samples at 4°C for 10 minutes at 1,000 x g. Remove serum and assay immediately or store samples on ice for up to 6 hours before assaying. Aliquots of serum may also be stored at below -70°C for extended periods of time. Avoid repeated freeze-thaw cycles.

Plasma: Collect plasma using EDTA-Na₂ as the anticoagulant. If possible, collect the plasma into a mixture of EDTA-Na₂ and Futhan5 to stabilize the sample against spontaneous *in vitro* complement activation. Immediately centrifuge samples at 4°C for 15 minutes at 1,000 x g. Assay immediately or store samples on ice for up to 6 hours before assaying. Aliquots of plasma may also be stored at below -70°C for extended periods of time. Avoid repeated freeze-thaw cycles.

Note: Heparin and Citrate plasma has not been validated for use in this assay.

Other biological samples: MBL has not tested.

(*e.g.* Remove any particulates by centrifugation and assay immediately or aliquot and store samples at below -70°C. Avoid repeated freeze-thaw cycles. Individual users should determine appropriate conditions when using other types of samples.)

CircuLex™



For Research Use Only, Not for use in diagnostic procedures

Detailed Protocol

The MBL Research Product **CircuLex S100A8/MRP8 ELISA Kit** is provided with removable strips of wells so the assay can be carried out on separate occasions using only the number of strips required for the particular determination. Since experimental conditions may vary, an aliquot of the Human S100A8 Standard within the kit, should be included in each assay as a calibrator. Disposable pipette tips and reagent troughs should be used for all liquid transfers to avoid cross-contamination of reagents or samples.

Preparation of Working Solutions

All reagents need to be brought to room temperature prior to the assay. Assay reagents are supplied ready-to-use, with the exception of 10X Wash Buffer, 20X HRP-conjugated Detection Antibody and Human S100A8 Standard.

- 1. Prepare a working solution of Wash Buffer by adding 100 mL of the **10X Wash Buffer** to 900 mL of deionized (distilled) water (ddH₂O). Mix well. Store at 4°C for two weeks or -20°C for long-term storage.
- 2. Prepare HRP conjugated Detection Antibody by <u>20-fold</u> diluting **20X HRP-conjugated Detection Antibody** with Conjugate Dilution Buffer at the time of assay*.

*Prepare appropriate volume for your assay. Discard any unused HRP-conjugated Detection Antibody after diluted.

 Reconstitute Human S100A8 Standard with X* μL of Dilution Buffer. The concentration of the human S100A8/MRP8 in vial should be <u>50 ng/mL</u>, which is referred as a Master Standard of human S100A8/MRP8.

*The amount is changed depending on lot. See the real "User's Manual" included in the kit box.

Prepare Standard Solutions as follows:

Use the **Master Standard** to produce a dilution series (below). Mix each tube thoroughly before the next transfer. The 5,000 pg/mL standard (Std.1) serves as the highest standard. The Dilution Buffer serves as the zero standard (Blank).

	Volume of Standard	Dilution Buffer	Concentration
Std.1	60 μL of Master Standard (50 ng/mL)	540 μL	5,000 pg/mL
Std.2	300 μL of Std. 1 (5,000 pg/mL)	300 µL	2,500 pg/mL
Std.3	300 μL of Std. 2 (2,500 pg/mL)	300 µL	1,250 pg/mL
Std.4	300 μL of Std. 3 (1,250 pg/mL)	300 µL	625 pg/mL
Std.5	300 μL of Std. 4 (625 pg/mL)	300 µL	313 pg/mL
Std.6	300 μL of Std. 5 (313 pg/mL)	300 µL	156 pg/mL
Std.7	300 μL of Std. 6 (156 pg/mL)	300 µL	78.1 pg/mL
Blank	-	300 µL	0 pg/mL

Note: Do not use a Repeating pipette. Change tips for every dilution. Wet tip with Dilution Buffer before dispensing. Unused portions of Master Standard should be aliquoted and stored at below -70°C immediately. Avoid multiple freeze and thaw cycles.



User's Manual



For Research Use Only, Not for use in diagnostic procedures

Sample Preparation

Dilute samples with **Dilution Buffer**.

• Serum and plasma samples may require a 100-fold dilution *e.g.* 3μ L of sample + 297 μ L of Dilution Buffer

Assay Procedure

- 1. Remove the appropriate number of microtiter wells from the foil pouch and place them into the well holder. Return any unused wells to the foil pouch, refold, seal with tape and store at 4°C.
- 2. Dilute samples with Dilution Buffer. (See "Sample Preparation" above.)
- 3. Pipette 100 μL of Standard Solutions (Std1-Std7, Blank) and diluted samples in duplicates, into the appropriate wells.
- 4. Incubate the plate <u>at room temperature (ca.25°C) for 60 minutes</u>, shaking at ca. 300 rpm on an <u>orbital microplate shaker</u>.
- 5. Wash 4-times by filling each well with Wash Buffer (350 μL) using a squirt bottle, multi-channel pipette, manifold dispenser or microplate washer.
- 6. Add 100 μL of HRP conjugated Detection Antibody into each well.
- 7. Incubate the plate <u>at room temperature (ca.25°C) for 60 minutes</u>, shaking at ca. 300 rpm on an <u>orbital microplate shaker</u>.
- 8. Wash 4-times by filling each well with Wash Buffer (350 μL) using a squirt bottle, multi-channel pipette, manifold dispenser or microplate washer.
- 9. Add 100 μL of Substrate Reagent. Avoid exposing the microtiter plate to direct sunlight. Covering the plate with e.g. aluminum foil is recommended. Return Substrate Reagent to 4°C immediately after the necessary volume is removed
- 10. Incubate the plate <u>at room temperature (ca. 25°C) for 10-20 minutes</u>, shaking at ca. 300 rpm on an <u>orbital microplate shaker</u>. The incubation time may be extended up to 30 minutes if the reaction temperature is below than 20°C.
- 11. Add 100 μ L of Stop Solution to each well in the same order as the previously added Substrate Reagent.
- 12. Measure absorbance in each well using a spectrophotometric microplate reader at dual wavelengths of 450/540 nm. Dual wavelengths of 450/550 or 450/595 nm can also be used. Read the microplate at 450 nm if only a single wavelength can be used. Wells must be read within 30 minutes of adding the Stop Solution.
- **Note-1:** Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Buffer by aspirating or decanting. Invert the plate and blot it against clean paper towels.



User's Manual



For Research Use Only, Not for use in diagnostic procedures

- **Note-2:** Reliable standard curves are obtained when either O.D. values do not exceed 0.2 units for the blank (zero concentrations), or 2.8 units for the highest standard concentration.
- **Note-3**: If the microplate reader is not capable of reading absorbance greater than the absorbance of the highest standard, perform a second reading at 405 nm. A new standard curve, constructed using the values measured at 405 nm, is used to determine human S100A8/MRP8 concentration of off-scale samples. The readings at 405 nm should not replace the on-scale readings at 450 nm.

Calculations

Average the duplicate readings for each standard, control, and sample and subtract the average zero standard optical density. Plot the optical density for the standards versus the concentration of the standards and draw the best curve. To determine the human S100A8/MRP8 concentration of each sample, first find the absorbance value on the y-axis and extend a horizontal line to the standard curve. At the point of intersection, extend a vertical line to the x-axis and read the corresponding human S100A8/MRP8 concentration. If the samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.

- 1. The dose-response curve of this assay fits best to a sigmoidal 4-parameter logistic equation. The results of unknown samples can be calculated with any computer program having a 4-parameter logistic function. It is important to make an appropriate mathematical adjustment to accommodate for the dilution factor.
- 2. Most microtiter plate readers perform automatic calculations of analyte concentration. The calibration curve is constructed by plotting the absorbance (Y) of calibrators versus log of the known concentration (X) of calibrators, using the four-parameter function. Alternatively, the logit log function can be used to linearize the calibration curve (i.e. logit of absorbance (Y) is plotted versus log of the known concentration (X) of calibrators).

Measurement Range

The measurement range is 78.1 pg/mL to 5,000 pg/mL. Any sample reading higher than the highest standard should be diluted with Dilution Buffer in higher dilution and re-assayed. Dilution factors need to be taken into consideration in calculating the human S100A8 concentration.

Troubleshooting

- 1. All samples and standards should be assayed in duplicate, using the protocol described in the **Detailed Protocol**. Incubation times or temperatures significantly different from those specified may give erroneous results.
- 2. Poor duplicates, accompanied by elevated values for wells containing no sample, indicate insufficient washing. If all instructions in the **Detailed Protocol** were followed accurately, such results indicate a need for washer maintenance.
- 3. Overall low signal may indicate that desiccation of the plate has occurred between the final wash and addition of Substrate Reagent. <u>Do not allow the plate to dry out</u>. Add Substrate Reagent immediately after wash.



User's Manual



For Research Use Only, Not for use in diagnostic procedures

Reagent Stability

All of the reagents included in the MBL Research Product **CircuLex S100A8/MRP8 ELISA Kit** have been tested for stability. Reagents should not be used beyond the stated expiration date. Upon receipt, kit reagents should be stored at 4°C, except the reconstituted Human S100A8 Standard must be stored at below -70°C. Coated assay plates should be stored in the original foil bag sealed by the zip lock and containing a desiccant pack.

Assay Characteristics


1. Sensitivity

Twenty-one assays were evaluated and the minimum detectable dose (MDD) of human S100A8/MRP8.

The MDD (defined as such a concentration of human S100A8/MRP8 giving absorbance higher than mean absorbance of blank* plus three standard deviations of the absorbance of blank: A blank + 3SD blank) is better than 43.4 pg/ml of sample.

* Dilution Buffer was pipetted into blank wells

Typical standard curve

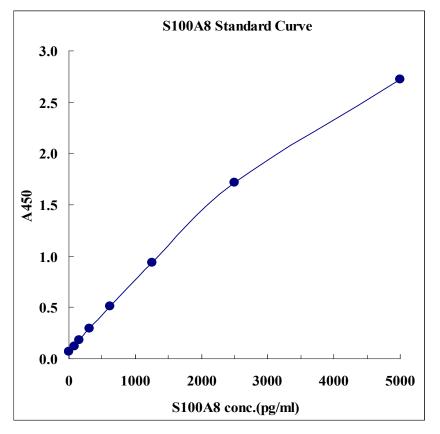






User's Manual For Research Use Only, Not for use in diagnostic procedures

2. Specificity

The antibodies in the human S100A8/MRP8 ELISA Kit are highly specific of human S100A8, with no detectable cross-reactivity to all other S100 proteins.

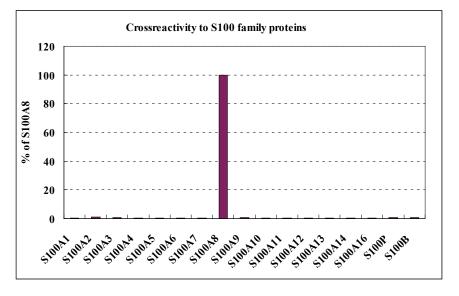

3. Precision

Intra-assay Precision (Precision within an assay)

Four samples of known concentration were tested sixteen times on one plate to assess intra-assay precision.

• Intra-assay (Within-Run, n=16) CV=3.6-4.1 %

Human S100A8/MRP8 conc. (ng/mL)

	Serum-1	Serum-2	Serum-3	Serum-4
1	121.4	144.4	196.5	348.8
2	120.8	142.8	189.2	328.6
3	125.5	148.3	199.9	344.3
4	124.5	151.7	202.7	351.7
5	121.8	151.2	194.8	345.1
6	120.3	147.3	199.5	353.3
7	117.5	152.0	186.7	337.0
8	113.7	133.5	174.9	320.1
9	118.6	142.4	203.2	366.5
10	117.9	146.9	195.7	328.4
11	117.2	148.5	200.3	349.8
12	120.7	152.8	201.1	334.3
13	120.5	153.0	199.9	331.9
14	120.7	147.9	193.2	341.0
15	113.1	141.9	192.9	320.6
16	108.4	143.2	194.4	315.3
max.	125.5	153.0	203.2	366.5
min.	108.4	133.5	174.9	315.3
mean	118.9	146.7	195.3	338.5
SD	4.3	5.2	7.2	14.0
CV(%)	3.6	3.6	3.7	4.1

Cat#: CY-8061

CircuLex[™]



User's Manual

For Research Use Only, Not for use in diagnostic procedures


Inter-assay Precision (Precision between assays)

Five samples of known concentration were tested in four separate assays to assess inter-assay precision.

	Serum-5	Serum-6	Serum-7	Serum-8	Serum-9
1	18.8	16.1	141.1	188.3	321.1
2	20.5	18.1	145.4	185.6	335.5
3	19.4	18.5	137.9	185.3	343.5
4	18.7	16.3	126.1	168.3	320.4
Max.	20.5	18.5	145.4	188.3	343.5
Min.	18.7	16.1	126.1	168.3	320.4
Mean	19.4	17.2	137.6	181.9	330.1
S.D.	0.8	1.2	8.3	9.2	11.3
C.V.(%)	4.3	7.1	6.0	5.0	3.4

4. Linearity

To assess the linearity of the assay, samples containing and/or spiked with high concentrations of human S100A8/MRP8 were serially diluted with the Dilution Buffer to produce samples with values within the dynamic range of the assay.

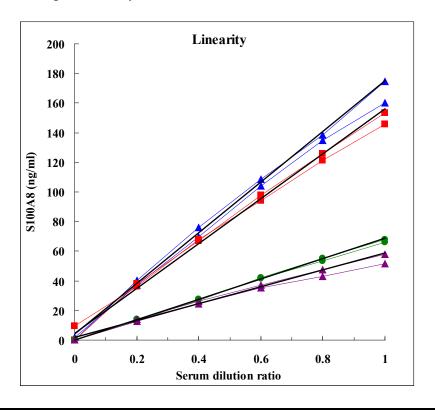






For Research Use Only, Not for use in diagnostic procedures

Example of Test Results

Fig.1 Concentrations of S100A8/MRP8 in Crohn's disease patients' sera (n=20) and healthy Japanese volunteers' sera (n=29)

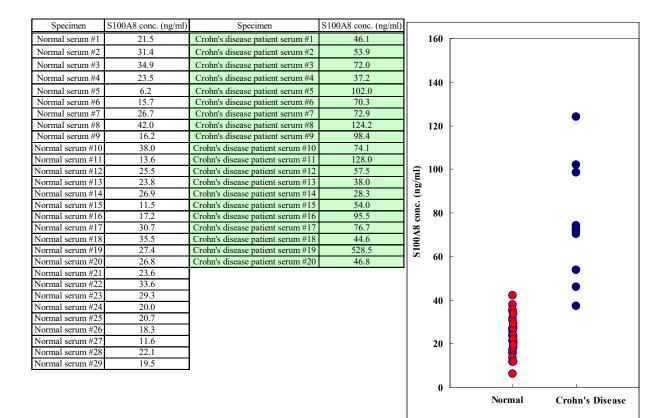



User's Manual



For Research Use Only, Not for use in diagnostic procedures

References

- 1. Odink, K. et al. (1987) *Nature* **330**, 80-82.
- 2. Kerkhoff, C., Klempt, M. & Sorg, C. (1998) Biochim. Biophys. Acta 1448, 200-211.
- 3. Sorg, C. (1992) Behring Inst. Mitt. 91, 126-137.
- 4. Foell, D., Hernandez-Rodriguez, J., Sanchez, M., Vogl, T., Cid, M. C., and Roth, J. (2004) *J. Pathol.* **204**, 311–316
- 5. Nacken, W., Roth, J., Sorg, C., and Kerkhoff, C. (2003) Microsc. Res. Tech. 60, 569-580
- 6. Foell, D., Frosch, M., Sorg, C., and Roth, J. (2004) Clin. Chim. Acta. 344, 37-51
- 7. Rammes, A., Roth, J., Goebeler, M., Klempt, M., Hartmann, M. & Sorg, C. (1997) *J. Biol. Chem.* 272, 9496–9502.
- 8. Xu, K., T. Yen, C. L. Geczy. (2001) J. Immunol. 166: 6358.
- 9. Grimbaldeston, M. A., C. L. Geczy, N. Tedla, J. J. Finlay-Jones, P. H. Hart. (2003) J. Invest. Dermatol. 121: 1168.
- Harrison, C. A., M. J. Raftery, J. Walsh, P. Alewood, S. E. Iismaa, S. Thliveris, C. L. Geczy. (1999) *J. Biol. Chem.* 274: 8561.

For more information, please visit our web site. https://ruo.mbl.co.jp/

MANUFACTURED BY

MEDICAL & BIOLOGICAL LABORATORIES CO., LTD. URL: https://ruo.mbl.co.jp E-mail: support@mbl.co.jp

CycLex/CircuLex products are supplied for research use only. CycLex/CircuLex products and components thereof may not be resold, modified for resale, or used to manufacture commercial products without prior written approval from MBL. To inquire about licensing for such commercial use, please contact us via email.